


3 Easy Steps

1

Connect CPAP tubing inside the bag. Put in mask & chamber*. Zip it closed.


Unzip the bag and find the round valve inside. Detach the tubing from your CPAP machine and attach it to this valve. Put all CPAP tubing, mask, and dry water chamber inside the bag, then zip the bag closed.

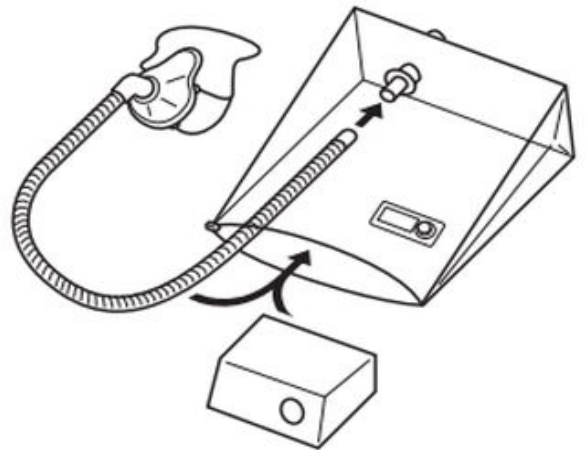

2

Connect Sleep8 outside the bag.

Slide your Sleep8's triangular port onto the triangular valve on the outside of the bag. It should be a snug fit. The device will only work once there's a secure connection.

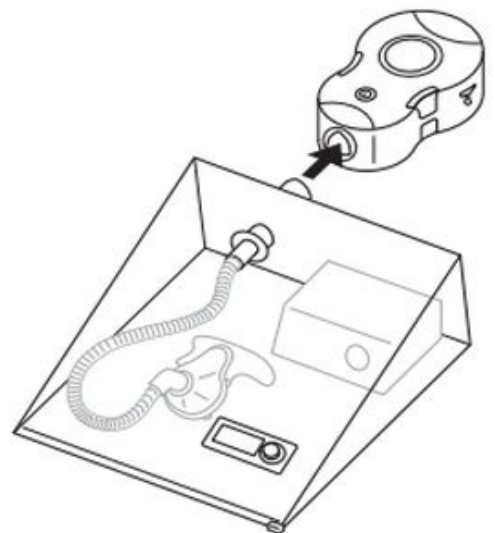

3

Press start & relax.

Press the power button until you hear a beep and see the Sleep8 screen light up. The screen will count down until the sanitization is complete. It's normal to hear the fan oscillate and pause for short periods. When done, your Sleep8 will automatically turn off. We recommend leaving the bag closed for 1 hour after the cleaning is completed. Then unzip the bag to find your sanitized CPAP ready to use.



Charging your Sleep8 device.

When you notice the battery is low, plug it in using the Sleep8 charger included. You cannot overcharge your Sleep8.

When the device is plugged in using the Sleep8 charger, the device will beep once, the screen will be dim, and the battery icon will show the charging progress. Once the icon shows 4 full bars, your Sleep8 is fully charged and ready for action.



Visit [MySleep8.com](https://mysleep8.com) for a video tutorial and other resources.